Product Data Sheet

ClimaCORR[®] CC 1000-TL (+80°C)

Relevant Test Standards:

- Cyclic Climate Tests:
- -DIN EN ISO 11997-1:2006 Cycle B (previously VDA 621-415) -VW PV 1210 (climate module
- required)
- -Nissan CCT I
- -Renault ECC1 (climate module required)
- -IEC-60068-2-52 KB

Salt Spray Test:

- -DIN EN ISO 9227 -DIN 50942, DIN 53167 -ASTM B 117-73, ASTM B 287-74 -ASTM B 368-68 -ISO 7253 ISO 3678 -BS 1224, BS 2011, BS3900 F4 -BS 3900 F12 -BS 5466 Part I, BS 5466 Parts 2 + 3
- -NFX 41002,
- -AS 21331 Section 3.1 -SIS 1841190
- -JIS Z 2371

Condensation water test:

-ISO 6270-2 -BS 3900F2, BS 3900 F15 -ASTM D2242



Figure 1 Patented Controlled Water Condensation (CWC) system

Order Information

Model: ClimaCORR[®] CC 1000-TL Article number: - V.745.672.470

Sales & Support: +49 5205 87963 0 Monday to Friday 8:00 am – 17:00 pm

VLM GmbH Heideblümchenweg 50

33689 Bielefeld

info@vlm-labtec.com www.vlm-labtec.com

Specification subject to changes Pictures might differ from original

1 CC 1000-TL | VLM GmbH



Applicable (Corporate) Cyclic Climate Tests

- DIN EN ISO 11997-1:2006 Cycle B (previously VDA 621-415)
- VW PV 1210 (external climate module required)
- Nissan CCT1
- Renault ECC1 (external climate module required)
- IEC-60068-2-52 KB

Product Description

Compact front loading model of the test chamber developed for the automated conduction of advanced cyclic corrosion tests pursuant to the most of the corporate standards of the international automotive industry i.e. DIN EN ISO 11997-1:2006 Cycle B (previously VDA 621-415), VW PV 1210, Nissan CCT1, Renault ECC1 as well as the standard Salt Spray test acc. to ISO 9227 and High Humidity 6270-2.

Customer Benefits

- Operating temperature range from ambient up to +80°C
- Advanced steel design of the chamber allows for fast temperature and humidity transient times making this test chamber perfectly suitable for accelerated corrosion tests – these result in time and energy savings during the operational lifetime
- The patented VLM technology allows the best possible reproducibility of the climate conditions being created in the test chamber regardless the environmental conditions and geographical location same type of VLM corrosion test chamber will produce same test results in any location in the world. In total 9 different sensors are monitoring all relevant test parameters being temperature, humidity, air flow during the drying phase and the purity of the demineralized water
- The test chamber made of steel is more robust and less susceptible for damages compared to the competitive products made of glass reinforced plastic
- Lower cost of ownership compared to the competitive products where the test chamber is made of glass reinforced plastic (shorter test periods, better energy efficiency, easier for service and maintenance, longer life cycle, more resistive to mechanical damages)
- User friendly, menu guided and intuitive control system based on a touch screen in which most common standard corrosion tests are already preconfigured

Product Data Sheet

ClimaCORR[®] CC 1000-TL (+80°C)

-Nissan CCT I

-BS 3900 F12


-NFX 41002.

-SIS 1841190

-JIS Z 2371

-ISO 6270-2

-ASTM D2242

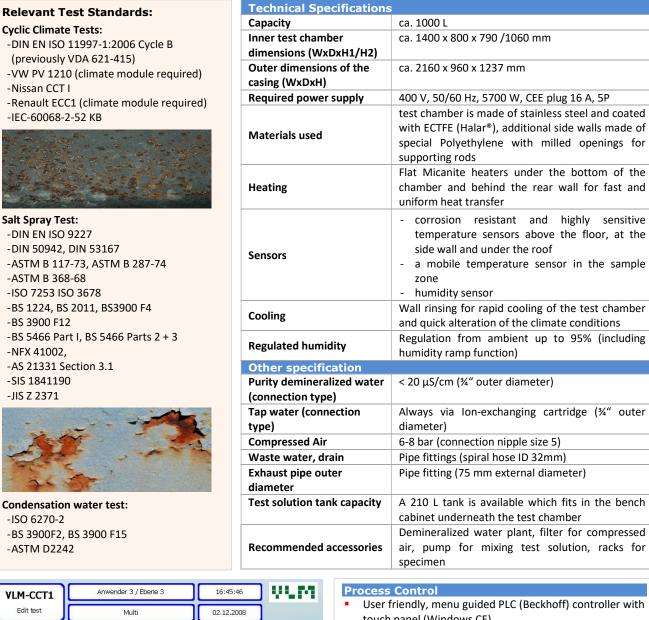

PR was	and a state of the state of the state of the state	ocken = Cha	mber dry CON Fill Warm = Fillin Steps	g Chamber	Segments
1	Salznebel 50	1	Salznebel 50°	1	Sprung
2	CON KK 80	2	CON Füllwarm 50°	2	Zyklusende
3	CON KEW AHT	3	Kondenswasser 50°	3	Salznebel
4	CON KTW AT	4	Warmluft / Belüften	4	Salznebel CASS
5	Prohesion	5	Sprung zu 1 / 11x	5	Salznebel Prohesion
6	Belüften	6	Zyklusende zu 1 / Ox	6	Kondenswasser
7	Warmluft	7	T	7	Belüften/PR Wasser 🚽
8	CCT 1	4	►	•	►.
9 Nai	ECC1 The [max. 20 Characters]	Veeks Days: Hours	0 🕹 Ĩ		Save Segment
	Warmluft	Minute	es: O		

Figure 2 Menu screen with test management functions



CC 1000-TL | VLM GmbH 2

- User friendly, menu guided PLC (Beckhoff) controller with touch panel (Windows CE)
- Ethernet interface for connecting to local intranet (features remote control and test data transfer)
- Restricted access for operators to three different operating levels
- Memory storage with the capacity to hold of 50 test programs, most common corrosion tests already preconfigured in the factory, users can easily program own tests
- All tests data and relevant test parameters are recorded and can be via Ethernet exported to a PC for further processing (e.g. Microsoft Excel)
- The Beckhoff PLC controls remotely the external climate module for the tests which require such module (typically PV 1210 and Renault ECC1)
- The Beckhoff PLC controls remotely the powering of test specimen (switch ON and OFF)
- Alarm table and history log
- Full overview of all digital and analog inputs / outputs

Version: v4/05.02.2015

Product Data Sheet

ClimaCORR[®] CC 1000-TL (+80°C)

Accessories included:

- 7 rods for supporting test specimen
- 2 m exhaust hose Ø 75 mm
- 2 m drain water hose Ø 32 mm
- 1 female connector for compressed air hose (size no. 5)
- 1 bag sodium chloride special quality for corrosion testing, accurate weight 5000g
- 1 set of measuring cylinders for measuring the fall our rate
- 1 spare filter for the test solution in the reservoir



Operating System Salt Spray (according to ISO 9227)

- Electronically controlled diaphragm pump for accurate flow control of the test solution with automatic monitoring of the presence of the air in the pump system
- Adjustable high performance acrylic spray nozzle for an optimal generation and distribution of the salt fog
- Option (standard for ClimaCORR[®] types): Second nozzle for spraying/washing of the test chamber with demineralized water after the salt spray test is finished
- Transparent humidifier made of heat resistant Borosilicate glass with a replaceable PE air filter for creating fine air bubbles in order to fully saturate the compressed air with humidity
- Automatic refilling of the test chamber for maintaining the constant water level in the trough
- Safety air pressure valve on the humidifier
- Over temperature protection of the humidifier
- Air purge function for rapid evacuation of the salt fog out of the chamber (after a salt spray test is finished) and prior to opening the door
- Accessories: 210 L PE storage reservoir for salt solution, level indicator is an option

Operating system Regulated Humidity

- The humidity in the test chamber at temperatures below 50°C is achieved by injecting the compressed air saturated with humidity
 into the test chamber through the spray nozzle
- The humidity in the test chamber at temperatures above 50°C is achieved with the special steam generator (patent pending) capable of injecting controllable quantities of steam into the test chamber
- Humidity is monitored by a high performance humidity sensor which is pneumatically inserted into the test chamber only when the operating system Controlled Humidity is active
- Operating range for controlled humidity (see Figure 1): from environmental up to 98% RH at 35°C
- Accuracy of the humidity regulation is ±3% over the full control range and ±1% in the high humidity control range (±90% RH)
- For humidity below environmental (typically <50% RH) an external climate module is required (e.g. to achieve 20% RH at 35°C as specified by ECC1 test)

Operating system High Humidity with Controlled Water Condensation - CWC (according to ISO 6270-2 CH)

- CWC system (Figure 1) is the patented VLM technology which regulates the temperature gradient of exactly ΔT=1°C between the bottom and the roof of the test chamber – this is essential for an optimal condensation process in the test chamber at 100% RH regardless the environmental conditions outside the test chamber
- Flat heaters under the bottom of the chamber and behind side walls allow for uniform and rapid heating of the chamber
- Temperature stability in the chamber ± 0,5°C
- Air fan with adjustable rotation speed for controllable drying of specimen in the Drying Phase; all ClimaCORR[®] types feature an
 extra air flow sensor which feeds the information about the air flow through the chamber back to the PLC controller which
 regulates the fan rotation speed according to common standard corrosion tests are already preconfigured

Operating system Forced Air Drying / Aeration

- Warm air up to +60°C for forced air drying phase
- Adjustable air distribution system for uniform drying of the test specimens
- Fan with controllable rotation speed and sensor for air flow (sensor standard only for ClimaCORR[®] types) for automatic regulation of the air flow
- Temperature control range: from ambient up to +60°C (accuracy ±1°C)
- Transition time from Dry 60°C to "Wet" Humidity 50°C > 95% (not condensing humidity) within 15 min
- Transition time from Wet 50°C to Spray 35°C within 20 min
- Transition time from Spray 35°C to Dry 60°C within 20 min
- Option: Air conditioning module to provide standard climate acc. to DIN 50014: ambient air at 23°C ±2°C and 50% ±5% RH
- Option: Air conditioning module to provide standard climate for conducting Renault ECC1 Test (30°C ±2°C and 20% ±5% RH)